



MEDIOSTAR

NEXT PRO FAMILY
HYBRID TECHNOLOGY

THE UNMATCHED STARS FOR
PROFESSIONAL HAIR REMOVAL



THE NEW

MEDIOSTAR NEX~~T~~ PRO FAMILY

BEYOND PROFESSIONAL HAIR REMOVAL

MEDIOSTAR HYBRID NEX~~T~~ PRO XL

The largest spot size available on the market and treatments faster than ever.

Try today the power of an extraordinary hybrid hair removal experience.

>>>

MEDIOSTAR NEX~~T~~ PRO ALEX

The efficiency of the Alexandrite wavelength meets the latest diode technology.

A unique combination that allows professional results even on thinner and lighter hairs. >>>

MEDIOSTAR HYBRID NEX~~T~~ PRO

The best choice for users who are starting in the hair removal business and want to get on board with the latest hybrid technology. Reliable, with no running costs, upgradable: simply ready for the success. >>>

MEDIOSTAR HYBRID LASER NEXT PRO XL

ALSO EQUIPPABLE WITH

+ STANDARD HANDPIECE

+ XL HANDPIECE

XTRAPOWERFUL
XTRAFAST
XTRALARGE



FASTER THAN EVER

The new MeDioStar NeXT PRO XL features the largest spot available on the market and offers an extraordinary hair removal experience at the highest speed. Thanks to **SmoothPulse Mode** the treatment is faster than ever and still painless. The new **360° skin cooling system** allows a highly efficient pre-cooling of the epidermis, protecting it from the risk of burns and allowing the treatment in any direction. This guarantees the best comfort for both the patient and the user. The **XL-S and XL-L tips** are specially designed for the treatment of both larger and smaller areas. The special hybrid **wavelength mix of 810 & 940 nm** allows to reach even the deeper follicles, for efficient and effective treatments.

ALL THE FLEXIBILITY THAT YOU NEED IN AN XL SIZE.

- + **LARGEST SPOT ON THE MARKET**
- + **TREAT ALL SKIN TYPES, ALL YEAR ROUND**
- + **HIGH SPEED TREATMENTS**
- + **360° SKIN COOLING SYSTEM INTEGRATED**
- + **NO CONSUMABLES**

* THE XL-TIPS ARE INTERCHANGEABLE

XL-S*

HANDPIECE

Spot size
3.0 cm²

Fluence
Up to 60 J/cm²

Frequency
Up to 12 Hz

XL-L*

HANDPIECE

Spot size
9.1 cm²

Fluence
Up to 20 J/cm²

Frequency
Up to 10.5 Hz



TREAT A BACK IN LESS THAN
FOUR MINUTES



MEDIOSTAR **NEXT PRO ALEX**

ALSO EQUIPPABLE WITH

+ **STANDARD HANDPIECE**

+ **XL HANDPIECE**



**ALEXANDRITE WAS GOOD.
WE MADE IT BETTER.**

THE NEXT GENERATION OF ALEXANDRITE LASERS

The MeDioStar NeXT PRO ALeX brings the **Alexandrite wavelength** and the **latest Diode Technology** together, combining the efficacy of the high melanin absorption of the Alexandrite with all the benefits of the Diode. Thanks to the highly selective 755nm wavelength, even thinner and lighter hairs can be treated with the best results. The **SmoothPulse Mode** allows fast and painless treatments and the **integrated skin cooling system** guarantees the highest safety. The proven reliability of the Diode Technology assures **the optimization of the treatment's procedure** (no time consuming heating up) and **no running costs**.

RECOVERING YOUR INVESTMENT HAS NEVER BEEN SO EASY.

- + HIGHLY EFFECTIVE ON THINNER AND LIGHTER HAIRS
- + TREAT UP TO SKIN TYPE V
- + INTEGRATED SKIN COOLING SYSTEM
- + SMOOTHPULSE MODE FOR PAINLESS TREATMENTS
- + FASTEST RETURN ON INVESTMENT

ALX HANDPIECE

Spot size
1.4 cm²

Fluence
Up to 35J/cm²

Frequency
Up to 12Hz



**ALEXANDRITE WAVELENGTH
WITH ALL THE BENEFITS OF
THE DIODE**



MEDIOSTAR HYBRID LASER NEXT PRO

ALSO EQUIPPABLE WITH
+ XL HANDPIECE

THE UPGRADABLE
SOLUTION.
PREPARE FOR
SUCCESS!



GET ON BOARD WITH THE LATEST TECHNOLOGY

The MeDioStar NeXT PRO is the most suitable solution for new users who want to start the hair removal business with the state-of-the-art technology. Thanks to the **highly reliable Hybrid Diode Technology** and the **user-friendly interface**, even less experienced users can become confident with the most advanced hair removal techniques in a few simple steps. The **SmoothPulse Mode** and the **integrated skin cooling system** allow fast, safe and effective treatments, with basically no pain. The special wavelength mix of 810 & 940 nm allow to reach even the deeper follicles, for efficient and effective treatments. On top of this, the device allows **multiple treatments** including **acne treatment**, **skin rejuvenation** and **skin whitening** and it is easily upgradable also for **vascular treatments**.

IT IS JUST EXACTLY WHAT YOU NEED FOR GETTING OFF TO A GOOD START.

- + **FAST, SAFE AND EFFECTIVE**
- + **TREAT ALL SKIN TYPES, ALL YEAR ROUND**
- + **INTEGRATED SKIN COOLING SYSTEM**
- + **SMOOTH PULSE MODE FOR PAINLESS TREATMENTS**
- + **MULTIPLE TREATMENT OPTIONS**

STANDARD HANDPIECE

Spot size
1.4 cm²

Fluence
Up to 47 J/cm²

Frequency
Up to 12 Hz



UPGRADABLE & FLEXIBLE:
THE BEST SOLUTION FOR YOUR BUSINESS



UPGRADES*

MAKE THE MOST OF YOUR DEVICE

BEYOND JUST HAIR REMOVAL VASCULAR HANDPIECE

MeDioStar NeXT PRO goes beyond just hair removal, with versatility and user-friendliness in other applications, such as the **treatment of leg veins**, hemangiomas and other ectatic vessels. With the special vascular handpiece featuring a wavelength of 940nm, it allows non-invasive and fast treatments, usually with no need for an anesthetic.



VASCULAR HANDPIECE

Spot size	0.12cm ²
Fluence	Up to 210J/cm ²
Frequency	Up to 4Hz

POWER IN YOUR HANDS HIGH POWER HANDPIECE

Patients with **darker skin types and very fine hairs** can be very difficult cases to deal with. The High Power handpiece offers the possibility to treat this condition with highest effectiveness and completely safely, thanks to the diode technology reaching up to 90J/cm² with very short pulses.



HIGH POWER HANDPIECE

Spot size	0.48cm ²
Fluence	Up to 90J/cm ²
Frequency	Up to 12Hz

THE PORTABLE SOLUTION TROLLEY

The specially designed trolley allows the easiest transportation and the **maximum space optimization**, including a storage space for footswitch and power cord in the back. A stylish additional feature, for those who like the maximum comfort.



*
OPTIONAL
AVAILABLE FOR
ALL DEVICES
OF THE FAMILY

ASCLEPION

COMPETENCE, EXPERIENCE, SUCCESS

Asclepion Laser Technologies has been operating for **more than 35 years** as a leader on the international medical laser scene as manufacturer of advanced laser systems for dermatology, aesthetic medicine and surgery. Based in the Jena Optical Valley, famous throughout the world as the homeland of the German optics industry, the company invests continuously in research and development of new technologies. Today, customers in **more than 60 countries** trust in the Made in Germany technology of Asclepion and its proven scientific expertise.

MeDioStar NeXT PRO, the 5th generation of Asclepion's hair removal technology, boasts an unmatched success all around the world, proof of its effectiveness, safety and high quality.



High quality means for us also a customer oriented philosophy. That's why we work hard every day in order to provide you not only the best technology but also a full range of services, to support you at 360°.

TRAINING & EDUCATION



Both regional and international trainings are permanently organized in our training center in Jena. That's how we provide our customers all the knowledge that they need for safely and effectively using our laser devices.

Find out more at:
www.asclepion.com/academy

MARKETING & WEBCLUB



A wide range of marketing tools is available on the Asclepion's WEBCLUB, an online platform designed for giving you a real time access to all the latest news and important documents for the success of your practice.

Find out more at:
www.asclepion.com/webclub

TECHNICAL SERVICE



Our highly specialized technical assistance, provided both directly and through a network of local distributors, follows the customer from the installation, to the assistance in case of defect, to the delivery of spare parts.

Find out more at:
www.asclepion.com/service


Asclepion

Laser Technologies

TECHNICAL SPECIFICATIONS

System Laser	MeDioStar Hybrid NeXT PRO High Power Diode Laser class 4	Spot size	9.1, 3.0, 1.4 or 0.48 cm ² (hair removal)* 0.12cm ² (vascular treatments)
Power	max. 2400W	Frequency	max. 12Hz*
Wavelength	800 – 950 nm*	Pulse length	8 – 400ms*
Fluence	max. 20J/cm ² (XL-L Handpiece) max. 60J/cm ² (XL-S Handpiece) max. 35J/cm ² (ALX Handpiece) max. 47J/cm ² (STANDARD Handpiece) max. 90J/cm ² (HIGH POWER Handpiece) max. 210J/cm ² (VASCULAR Handpiece)	Display	8.4" LCD
		Sizes	Table Top: 36.5x56x46.5cm With Trolley: 36.5x56x93cm
		Weight	ca. 35 kg

*According to handpiece. All specifications are subject to change without notice.

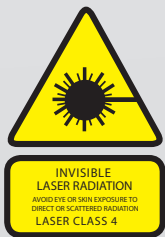

MADE IN GERMANY



PATENTS
DE 197 10 676
US 6,162,218



Management
System
EN ISO
13485:2012
www.tuv.com
ID 0000044235



Aesthetic - Medical - Foreign Trade Ltd/Co

MEDLASER AESTHETIC-MEDICAL-FOREIGN TRADE Ltd/Co

ANKARA : Gökkuşuğu Mahallesi 1204. Sokak No: 15 Cevizlidere - ÇANKAYA / TURKEY
Tel : +90 312 285 40 40 (Pbx)
Fax : +90 312 285 99 22
İSTANBUL : İttri Dede Sokak No: 16/1 Kızıltoprak / TURKEY
Tel : +90 216 349 33 49 (Pbx)
Web : www.medlaser.com.tr

Asclepion Laser Technologies GmbH | Brüsseler Str. 10 | 07747 Jena | Germany | www.asclepion.com

MEDIOSTAR NEXT PRO FAMILY

FOLLOW US ON

